
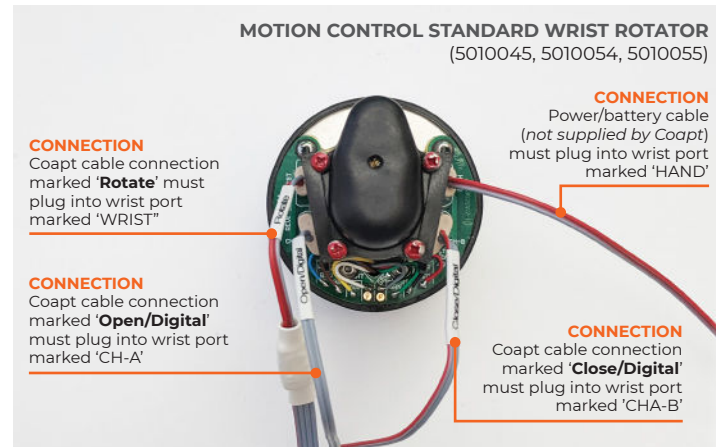


13

CLINICAL CONNECTIONS & SETTINGS GUIDE

I3 controls pronation/supination of powered wrists and direct pattern recognition command of multifunction grips for i-Limb Quantum* or VINCENTevolution3. (*with serial numbers M9618 or newer)

PART 1: PHYSICAL CONNECTIONS

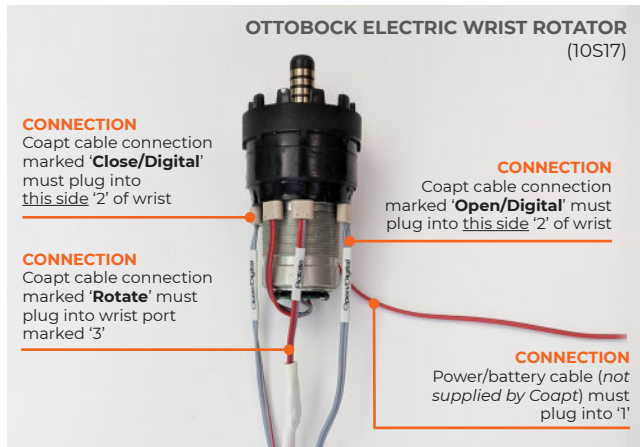

PART 2: SOFTWARE SETTINGS

WRISTS

OTTOBOCK ELECTRIC WRIST ROTATOR (NO MYOROTRONIC)

No settings required

MOTION CONTROL STANDARD WRIST ROTATOR (NO PROCONTROL)

No settings required

HANDS

I-LIMB QUANTUM*

Ensure hand settings are the factory defaults.

For proper operation of Gesture Control via Coapt, ensure 'CC' is selected as the i-mo trigger in the i-Limb app

VINCENTEVOLUTION3

Ensure hand settings are the factory defaults