


# 14

## CLINICAL CONNECTIONS & SETTINGS GUIDE

14 controls pronation/supination of Motion Control “6-band” powered wrist and direct pattern recognition command of multifunction grips for past-generation bebionic\*, TASKA, or Psyonic Ability Hand.

(\*with serial numbers BBHxxxQD, hand will require digital unlock code available upon request from Ottobock)

### PART 1: PHYSICAL CONNECTIONS

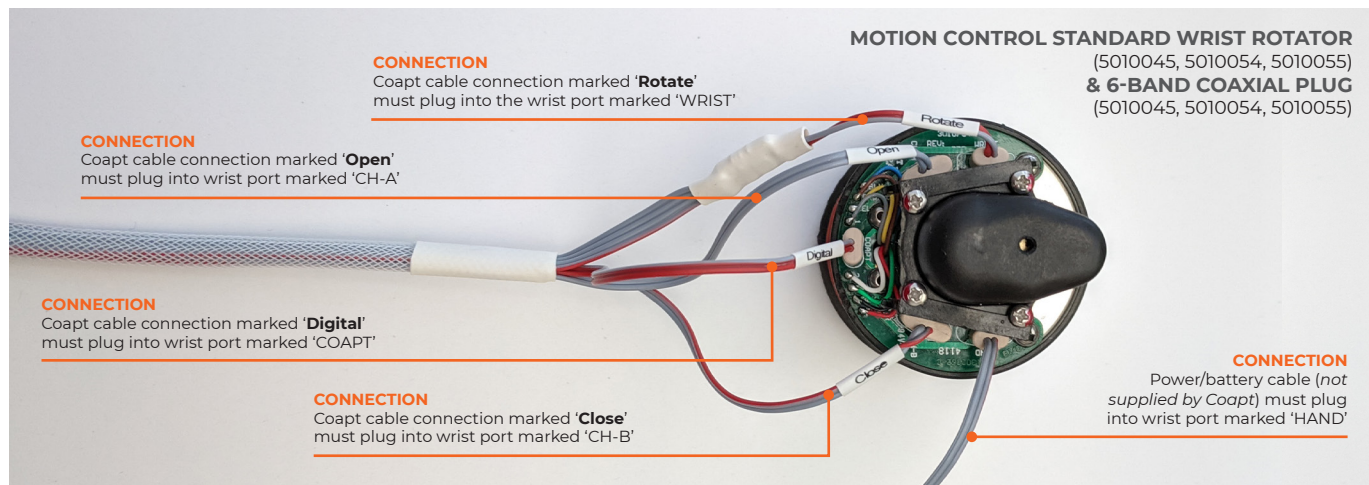

### PART 2: SOFTWARE SETTINGS

#### WRISTS

##### MOTION CONTROL STANDARD WRIST ROTATOR & 6-BAND COAXIAL PLUG

No settings required

#### HANDS

##### OTTOBOCK BEBIONIC\*

Control Strategy: Dual Site

Open / Close Strategy: Electrodes: A,B

Control Options: Electrodes

Control Response: Proportional

In DEVICE INFORMATION:

- Verify Firmware Rev: v304.032 or greater
- Click License Key
- Note the Serial Number, BT MAC and Contact Ottobock to obtain a License Key
- Click Enable once License Key is entered, then Close
- Press “Send to Hand” to confirm settings changes

In ELECTRODE SCREEN:

- Channel 1 “Blue” threshold: 10%, Channel 1 “Red” threshold: 90%, Channel 2 “Blue” threshold: 10%, Channel 2 “Red” threshold: 90%

Remember to press “Send to Hand” after making any changes

##### TASKA

TASKA settings will automatically be set upon each power-on

##### PSYONIC ABILITY HAND

Ensure hand settings are the factory defaults